


Belimed MST-V

Professional steam sterilisation of medical devices for CSSDs.





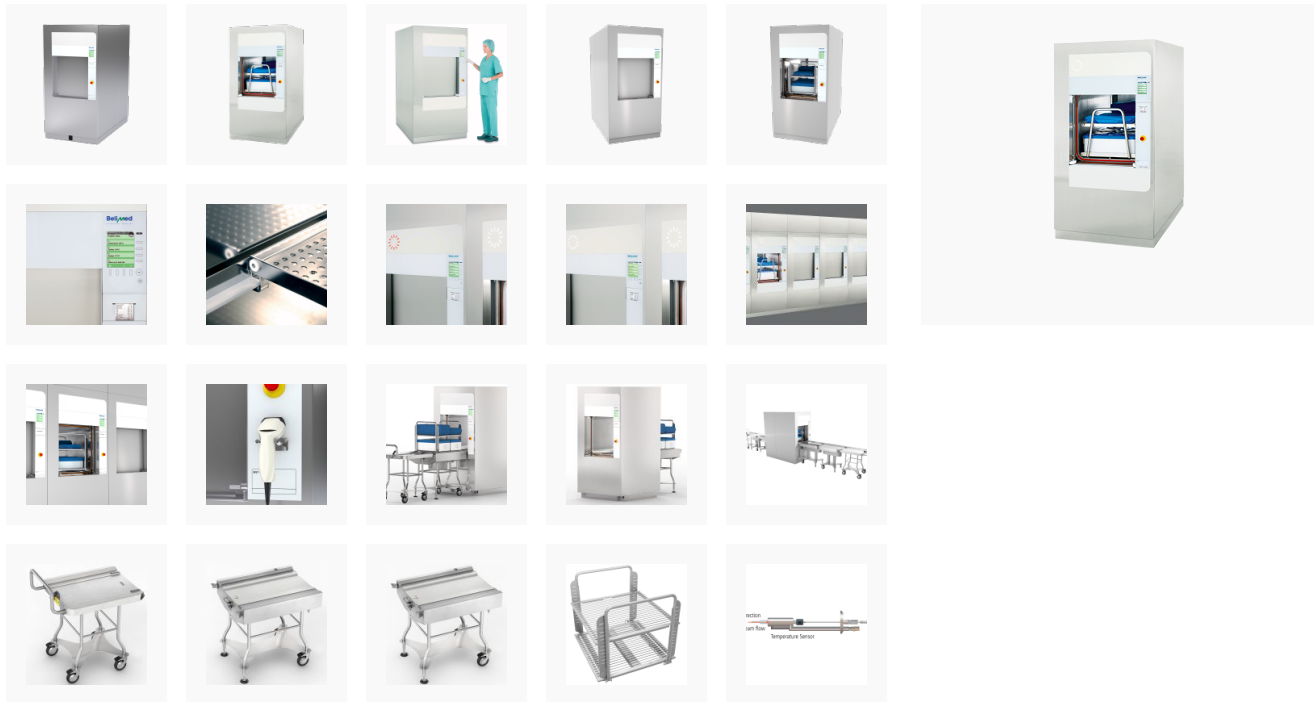









ARTG 197891

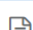

Features

- Belimed MST-V steam steriliser offers a large capacity throughput with a space-saving, efficient and ergonomic design, ideal for ensuring reliable and qualified sterilisation of medical devices in CSSD.
- Belimed's MST-V medical sterilisers user centric design combines excellent ergonomics with a loading height of 78 cm, lowest in this class with the touch-sensitive control panel also located at an ideal, ergonomic height.
- Clear, easy-to-see system process data can be done easily at a glance, from across the room with the Belimed MST-V steam steriliser.
- Vertical sliding doors and seamless surface made of solid stainless steel meets the highest standards in terms of hygiene.
- Belimed's MST-V control panel is located behind a glass plate, completely protecting it from soiling, making it quick and easy to clean.
- Belimed MST-V steam steriliser achieves maximum economy in terms of resource consumption. Innovative technology reduces water consumption, allowing for a higher inlet temperature of the cooling water and reduction of consumption is achieved by the water-saving system thanks to the integrated, extremely quiet liquid ring vacuum pump.
- Belimed MST-V is now available with optional integrated Bowie-Dick test resulting in no additional extra effort or materials required; it's electronic, so no consumables required, saving costs and is automated. The test can be programmed to preheat the chamber and run the test before the start of the CSSD shift.
- Belimed MST-V steam steriliser is perfect for maximum performance in small spaces. Offering capacity of 4 to 12 sterilisation units, its light weight and is only 99 cm wide, making it the most efficient and compact steriliser in its class.
- Belimed MST-V medical steriliser features a solid silicone door seal, with a milled solid stainless steel frame made of grade 1.4404 chrome-nickel steel, providing a much longer service life. Maintenance on the seal is fast and easy with Belimed providing a generous warranty on the door seal, backing its quality.

- Thanks to a wide range of versatile accessories, the degree of automation of the Belimed MST-V steriliser can be adjusted to your individual needs, thus enabling optimal process work flow.
- Belimed MST-V steam sterilisers offer easy front access to the technical equipment for quick maintenance. Front panelling components can be opened over the entire surface and the electric switch cabinet can be folded to the side, all assisting to save time, space and costs.
- Increase workflow efficiency, reliability and ease of use through Belimed's comprehensive and integrated automation systems for MST-V steam sterilisers.
- Belimed sterilizers are all engineered in Switzerland and are subjected to strict quality inspections to meet all the current standards and medical devices directives. That way they ensure maximum safety and quality and back the MST-V steriliser with generous warranty on the chamber and other components.
- Belimed SmartHub digital platform collects data from your machines and delivers them to an Instrument Tracking System like CensiTrac, offering simple transparency and provides you with the most important details of your Belimed machines at a glance.

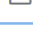Product Code	Product Description	Pack Size	Quote
107-7740	Belimed MST-V, 6-6-6, 1 Door	ea	
107-7761	Belimed MST-V, 6-6-9, 1 Door	ea	
107-7762	Belimed MST-V, 6-6-12, 1 Door	ea	
107-7763	Belimed MST-V, 6-6-6, 2 Door	ea	
107-7764	Belimed MST-V, 6-6-9, 2 Door	ea	
107-7765	Belimed MST-V, 6-6-12, 2 Door	ea	
107-7766	Belimed MST-V, 6-6-18, 2 Door	ea	
105-3678	Batch Cart to suit MST-V 6-6-6, without shelf	ea	
105-3679	Batch Cart to suit MST-V 6-6-9, without shelf	ea	
105-3680	Batch Cart to suit MST-V 6-6-12, without shelf	ea	
105-3711	Shelf for batch cart 6-6-6	ea	
105-3712	Shelf for batch cart 6-6-9	ea	
105-3713	Shelf for batch cart 6-6-12	ea	
105-3681	Transport Trolley to suit MST-V 6-6-6	ea	
105-3682	Transport Trolley to suit MST-V 6-6-9	ea	
105-3683	Transport Trolley to suit MST-V 6-6-12	ea	
TA16	NEW TS Automation System to suit MST-V Steam Sterilisers	ea	
1040274	Belimed Smart Hub	ea	

Technical Data

Dimensions (mm)	MST-V 6-6-6	MST-V 6-6-9
Product Outer (H x W x D) - 1 Door	1970 x 990 x 1095	1970 x 990 x 1395
Product Outer (H x W x D) - 2 Door	1970 x 990 x 1120	1970 x 990 x 1420
Product Chamber (H x W x D) 1 & 2 Door	660 x 660 x 700	660 x 660 x 1000
Chamber Volume (Sterilisation Units 330 x 300 x 600 mm)	4	6

Dimensions (mm)	MST-V 6-6-12	MST-V 6-6-18
Product Outer (H x W x D) - 1 Door	1970 x 990 x 1750	1970 x 990 x 2395
Product Outer (H x W x D) - 2 Door	1970 x 990 x 1775	1970 x 990 x 2420
Product Chamber (H x W x D) 1 & 2 Door	660 x 660 x 1350	660 x 660 x 2000
Chamber Volume (Sterilisation Units 330 x 300 x 600 mm)	8	12

Relevant Product

- [Belimed Elmasteam 8 Med](#)
- [Maxistream Steam Station](#)
- [Workstation Steam Box](#)
- [Belseal 95V](#)
- [Belseal 85V](#)
- [Belseal Pouchmaker](#)
- [Smeg WD2145](#)
- [Smeg WD1160](#)
- [Smeg WD4060](#)
- [Smeg WD4190](#)
- [UW565 Ultrasonic Cleaner](#)
- [Bandelin Sonorex](#)
- [Bandelin Sonomic](#)
- [Bandelin Triston](#)
- [Bandelin Sonoboard](#)

Applications

- CSD
- Healthcare Facilities
- Life Science

Resources

- Video
- video
- Video
- Video
- Video
- Brochure