

Soluscope Sprint AER

Introducing our state-of-the-art Automated Endoscope Reprocessor, designed to revolutionize the way you handle [endoscope reprocessing](#). This cutting-edge device boasts a fast cycle time, ensuring that your endoscopes are ready for use in record time. This innovation improves workflow efficiency in busy medical environments and helps to manage a high-volume hospital or a specialised day surgery/clinic with confidence. Let our reprocessor meet your needs with unparalleled speed and reliability.

Additionally, Soluscope Sprint AER features an enhanced pre-cleaning function, which meticulously prepares the endoscopes for the main cleaning cycle. The advanced pre-cleaning process effectively removes any residual debris and contaminants, paving the way for a thorough and efficient main cleaning cycle. The result is optimal cleanliness and hygiene, which is crucial for patient safety and infection control.


Our Automated Endoscope Reprocessor is fully compliant with industry standards, including the [GENCA](#) guidelines and [GESA](#) recommendations for endoscope reprocessing. These standards guarantee that our device meets the highest requirements for safety and performance, giving you peace of mind that your endoscopes are processed to the utmost standards of cleanliness and safety. Compliance with these guidelines means that our reprocessor not only meets but exceeds the stringent demands of modern healthcare environments, ensuring that your facility remains at the forefront of medical technology and patient care."

Experience the [future of endoscope reprocessing](#) with our innovative solution, designed to meet the rigorous demands of modern healthcare environments. With its fast cycle time and enhanced pre-cleaning function, our Automated Endoscope Reprocessor is the ideal choice for healthcare professionals seeking efficiency, reliability, and superior performance.

Invest in the best for your facility and ensure that your endoscopes are always ready for use, processed quickly and safely to the highest standards.

Take advantage of our [pay-as-you-go](#) option to make this advanced technology accessible and cost-effective for your facility.

ARTG 153397

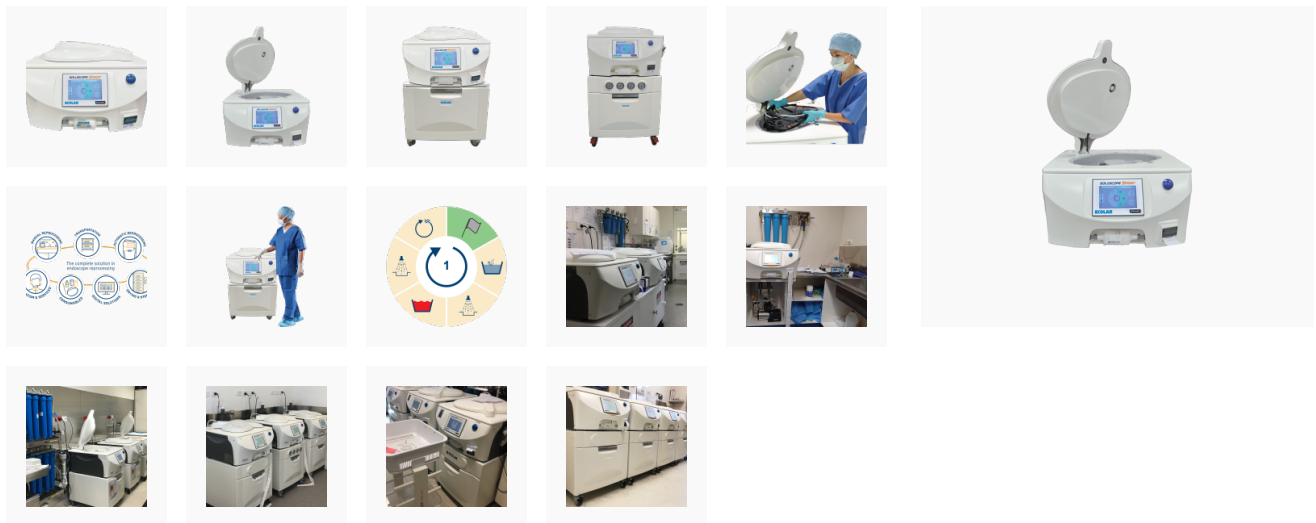


Features

Soluscope Sprint: Market Leader in Automated Endoscope Reprocessors (AERs)

- Soluscope Sprint is a market leader in flexible Automated Endoscope Reprocessors (AERs), backed by 25 years of endoscopy expertise. Our AERs comply with Australian Standards AS 5369:2023 and are registered with the TGA.
- Soluscope Sprint offers fast cycles and a complete clean disinfection cycle, preventing backlog of dirty scopes and enhancing operational efficiency.
- Soluscope Sprint features a 7" easy-to-use intuitive graphical touchscreen, enabling users to access all cycle and machine parameters quickly and easily, improving workflow.
- Soluscope Sprint eliminates manual cleaning procedures and alcohol injection flushes with our optional Enhanced Pre-

Cleaning function, ensuring thorough endoscope reprocessing.

- Soluscope Sprint is compatible with leading flexible endoscope brands and types, ensuring patient safety in Healthcare and improved infection control.
- Soluscope Sprint includes individual channel monitoring with alarms, ensuring all endoscope channels are fully cleaned and disinfected or sterilised during the reprocessing cycle, complying with GENCA and MTANZ guidelines.
- Soluscope Sprint features a full leak testing cycle with alarm, preventing fluid ingress and ensuring endoscope integrity.
- Soluscope Sprint offers an automated self-disinfection cycle using a second biocide, enhancing infection prevention.
- Soluscope Sprint provides the convenience of both high-level disinfection and sterilisation cycles in one machine, optimising healthcare operations.
- Soluscope Sprint’s barcode tracking system incorporates the tracking of the endoscope, operator, doctor, and patient, ensuring traceability and documentation.
- Soluscope Sprint offers the ability to reprocess EUS (ultrasound) endoscopes, expanding its versatility in medical device reprocessing.

| Product Code | Product Description | Pack Size | Quote |
|--------------|----------------------|-----------|---|
| SL-V3-PA-G | Soluscope Sprint AER | ea |  |

Relevant Product

- [Soluscope Sprint PT](#)
- [Soluscope DSC8000 Drying Cabinets](#)
- [Soluscope EP-Digital](#)
- [Censis ScopeTrac](#)
- [Asepti Epizyme Ultra, 2 L & 5 L](#)
- [Water Testing, Table 7.3](#)
- [Soluscope Service Contract](#)
- [PQ / Validation](#)
- [Smartline Sinks](#)
- [Smartline Benches](#)
- [Cleaning Consumables](#)
- [Procedure Consumables](#)

Applications

- Endoscopy
- Healthcare Facilities

Resources

- Brochure